Silk'n™

Beautiful Technology

Motion

H3220

- EN INSTRUCTIONS FOR USE
- DE BEDIENUNGSANLEITUNG
- FR INSTRUCTIONS D'UTILISATION
- NL GEBRUIKSIINSTRUCTIES
- ES INSTRUCCIONES DE USO
- PT INSTRUÇÕES DE UTILIZAÇÃO
- GR ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ
- TR KULLANIM KILAVUZU
- IT MANUALE OPERATIVO
- NO BRUKSANVISNING
- SV BRUKSANVISNING
- DA BRUGERVEJLEDNING

Legal Notice

Home Skinovations Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

SILK'N and the SILK'N logo are registered trademarks of Home Skinovations Ltd, Tavor Building, Shaar Yokneam POB 533, Yokneam 2069206, ISRAEL www.silkn.eu or www.silkn.com, info@silkn.com.

Table of contents

EN

1. Get to know your Silk'n Motion	3
2. Safety with Silk'n Motion	4
3. Contraindications & warnings	4
4. Do not use Silk'n Motion if:	6
5. Consult your physician before use if you have:	7
6. Safety tips.....	7
7. Possible side effects.....	8
8. Long-term hair removal the Silk'n way	10
9. Get started!.....	14
10. After treatment	16
11. Maintenance.....	16
12. Troubleshooting	17
13. Customer service.....	18
14. Warranty	18
15. Disposal	18
16. Specifications	19
17. Signs and symbols	20

1. Get to know your Silk'n Motion

Silk'n Motion is a light-based device for permanent hair reduction designed for home-use. It is a powerful electrical device and should be used with special attention to safety. Be sure to read all warnings, contraindications and safety precautions before use, and follow them strictly. We recommend you re-familiarize yourself with the instructions before each treatment.

1.1 Package contents

- Hair removal device
- Adapter with cable
- Storage case
- Cleaning cloth
- Warranty card
- User manual

1.2 Device description

Refer to the diagram at the front of this manual to see the location of each feature of your Silk'n Motion .

1. Energy level indicator lights (1-5)
2. Skin colour sensor light
3. Ready light
4. Control button
5. Skin colour sensor
6. Treatment surface
7. Adapter socket

Intended use

Silk'n Motion is intended to remove unwanted hair and effect long-term, or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs re-growing after a treatment regime. You may use Silk'n Motion to remove unwanted body hair. Ideal body areas include the underarms, bikini line, arms, legs, face (below the cheekbones), back, shoulders and chest.

Do not use this device for any purpose other than as described in this manual.

2. Safety with Silk'n Motion

With Silk'n Motion safety comes first

Home Pulsed Light™ (HPL) technology in Silk'n Motion provides superior safety with lower energy levels. HPL achieves long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The low energy reduces the risk of harm or complications, and contributes to your overall safety.

3. Contraindications & warnings

Silk'n Motion is not designed for everyone.

Hair removal by lasers or intense pulsed light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage when used on the face and neck. For more information visit www.silkn.eu or your local Silk'n domain.

 This device must not be used by the following people: children, adults with limited physical, sensory or psychological capacities, or those lacking experience or knowledge in how to use the device. Children should be supervised to ensure that they do not play with the device.

 Do not expose treated areas to the sun, and do not use energy levels 2-5 on tanned skin or after recent sun exposure (including artificial tanning machines). Refer to section "6.1 Exposure to the sun".

- ❗ Use below the cheekbones and not around the eyes.
- ❗ Do not use on areas with injected fills or toxins.
- ❗ Do not treat on tattoos or permanent makeup, dark brown or black spots, (such as large freckles, moles, or warts), nipples, genitals, or lips.
- ❗ Do not treat on active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- ❗ Abnormal skin conditions caused by diabetes, systemic or other metabolic diseases.
- ❗ Do not treat on an area that has a history of herpes outbreaks.
- ❗ Do not use if you are pregnant or nursing.
- ❗ Do not treat on areas where you may want hair to grow. The results are permanent.
- ❗ Do not treat the same area of skin more than once per hair removal session.
- ❗ Stop use immediately if your skin is very red, or if your skin blisters or burns.
- ❗ Use a lower energy level if you experience any discomfort.
- ❗ Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol), or acetone to clean the skin before use.

Keep your device away from water or liquids! Do not use near bathtubs, showers, basins or other vessels containing water. Do not place or store the device where it can fall, be pushed or placed anywhere near water. There is a risk of electrocution. Unplug immediately if the device comes into contact with any fluids, and contact Customer Service.

Remove all jewellery before use.

Never leave the device unattended when plugged into an outlet. Keep the AC/DC adapter accessible at all times, so it can be removed from the socket easily in case of emergency.

Do not attempt to modify, open or repair your device. You may expose yourself to dangerous electrical components or pulsed light energy. You may cause serious harm to yourself, and you will void your warranty.

4. Do not use Silk'n Motion if:

- The plug, attachments or accessories for this device are not recommended by Home Skinovations Ltd.
- The device has a damaged cord or plug. Keep the power cord away from heated surfaces.
- You see or smell smoke when it is in use.
- The *skin colour sensor*, *treatment surface*, or *device casing* is cracked or coming apart.
- The device is not working or appears to be damaged in any other way.

In such cases, contact Silk'n Customer Service.

5. Consult your physician before use if you have:

- A history of keloidal scar formation.
- Epilepsy.
- A disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- A history of skin cancer or areas of potential skin malignancies.
- Received radiation therapy or chemotherapy treatments within the past 3 months.
- Used physician-prescribed drugs on your skin within the past 6 months.
- Any other condition which, in your physician's opinion, would make it unsafe for you to be treated.

6. Safety tips

This section has information to help you use your device safely.

6.1 Exposure to the sun

- The lowest energy level (1) is designed for safe use even after sun exposure.
- If you use energy levels 2-5, avoid further exposure to the sun for 4 weeks before and 2 weeks after your treatment. This is because higher quantities of melanin expose the skin to greater risk when using any method of light-based hair removal. This applies to all skin types, even those that don't seem to tan quickly.
- We don't advise it, but if you have to go out in the sun after a treatment, be sure to apply sunscreen SPF30 or higher to the treated area. Do this for 2 weeks after a treatment.

6.2 Exposure of pulse light to the eyes

- Although you can use the device to treat facial hair, take extra care to avoid contact with the eyes, and use only below the cheekbones.

6.3 Wounds or infections

- Wait for the affected area to heal before further use.

7. Possible side effects

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. However every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore, it is important that you understand and accept the risks and complications that may occur with pulsed light hair removal systems designed for home use.

Minor skin discomfort

Home pulsed light hair removal is generally well-tolerated. However, some users feel mild discomfort during use. This is best described as a mild stinging sensation on the treated skin areas. The stinging sensation may last during or after treatment.

This means that you should use Silk'n Motion at a lower energy level, and if that doesn't help, you should stop using the device.

Skin redness

Your skin may become red immediately after using the device or within 24 hours. Redness generally clears up within 24 hours.

See your doctor if redness does not go away within 2 to 3 days.

Make sure you use the energy level that is comfortable for you and suitable for your skin colour.

Increased sensitivity of the skin

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Pigment changes

Silk'n Motion targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless there is a risk of temporary hyperpigmentation (increased pigment or brown discolouration) or hypo-pigmentation (whitening) to the surrounding skin. Risk of changes in skin pigmentation is higher for people with darker skin tones or heavily tanned skin. Usually discolouration or changes to skin pigment are temporary, and permanent hyperpigmentation or hypo-pigmentation are rare.

Excessive redness and swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but avoid exposure to the sun.

Blisters or burns

Cool the affected area and apply anti-burn cream. Wait until the area has healed completely before another treatment, and use a lower energy level.

If the skin does not return to normal within 3 days, consult your physician.

8. Long-term hair removal the Silk'n way

The efficiency of hair removal varies from person to person according to body area, hair colour, and how Silk'n Motion is used. A typical full hair growth cycle may take 18-24 months. During this time, multiple treatment sessions are required to achieve permanent hair removal.

A typical Silk'n Motion hair removal plan during a full hair growth cycle includes:

- Treatments 1-4: Plan two weeks apart.
- Treatments 5-7: Plan four weeks apart.
- Treatments 8+: Treat as needed, until desired results are achieved.

8.1 Expected treatment outcomes

Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

Skin Type ¹	Body Hair ²	Light Fluence [J/cm ²]	Part of the Body	Average Treatment Time	Effectiveness of Hair Removal
1-5	Dark blonde to brown	2.5-4	Lower or upper legs	Session 10-12, approx. 4 weeks apart.	60%
			Arms		50%
			Bikini line		60%
			Underarms		55%

Skin Type ¹	Body Hair ²	Light Fluence [J/cm ²]	Part of the Body	Average Treatment Time	Effectiveness of Hair Removal
1-5	Brown to black	2.5-4	Lower or upper legs	Sessions 8-10 approx. 4 weeks apart.	70%
			Arms		60%
			Bikini line		70%
			Underarms		65%

¹ Skin types are based on the Fitzpatrick scale of skin colours. Refer to the skin colour chart in the diagram at the front of the manual.

² Not effective on white body hair.

8.2 Energy Levels


The energy level determines the intensity of the light pulse delivered to your skin. As the energy level increases, so do hair removal results, as well as the risk of possible side effects and complications.

- The device has 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 5 *energy level indicator lights*.
- When the device is turned on, the energy level is set to level 1 automatically. Only one energy level indicator light will be on.
- To change the energy level, press the *control button* repeatedly to increase the energy level from 1 to 5. Once at level 5, press the *control button* again to reduce the energy level from 5 back to 1.

To determine the right energy level for your treatment area, refer to section "9.1 Before your first treatment".

8.3 Pulsing and gliding

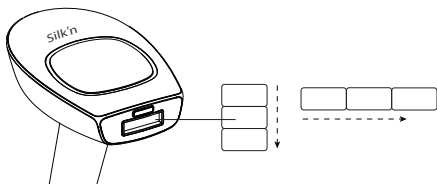Pulsing is the best technique for working at higher energy levels.

- Press the *treatment surface* against your skin. Work in rows (either across or down as shown in the image). Start at one end of each row and progress to the other end.
- The *treatment surface* creates temporary pressure marks on the treated area so you can see areas you have treated.
- After each pulse, move the *treatment surface* to the next spot. Keep doing this pulsing action for the area you are treating.
- Avoid gaps or overlaps as you use the device.

Gliding is the best technique for working at lower energy levels.

- Press the *treatment surface* firmly against your skin and slowly glide over your skin (either across or down as shown in the image). Make sure that the *treatment surface* is always touching your skin.
- At lower energy levels, the device takes only 1-1.5 seconds to recharge between pulses. This generates a nonstop sequence of flashes.
- Glide at a speed so that you avoid gaps or overlaps.
- ❗ Do not treat the same area of skin more than once per hair removal session.
- ❗ Stop use immediately if you have any discomfort.

8.4 Device sensors

Your Silk'n Motion has three sensors: for touch, movement, and colour.

- **Touch and Movement:** The device will only emit repeated pulses when it is in firm contact with your skin, and you keep moving it across your skin.
- **Colour:** Light-based hair removal on darker skin tones may result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypo-pigmentation). Your Silk'n Motion has a *skin colour sensor* that detects the colour of your skin automatically. The *skin colour sensor* checks your skin colour at the start of, and occasionally during, each treatment.

If Silk'n Motion detects that your skin tone is dark, and the energy level is set to 3 or higher, the device will auto-adjust and set the energy level back to 2.

8.5 What to expect when treating with Silk'n

For many people, using Silk'n Motion may be their first experience with a light-based device. During a treatment session it is normal to experience or feel:

- **Noise:** The cooling fan makes noise similar to a hair dryer. This is normal.
- **A popping sound:** With each flash of HPL™ light.
- **A sensation of warmth and tingling:** During each light pulse.
- **Some redness:** It is not uncommon to see some very mild, pink-like skin colouration during or just after a treatment. This is usually most noticeable around the hairs themselves.

Stop use immediately if you see your skin is very red.

Contact your physician immediately, if you experience any complication.

9. Get started!

To use your Silk'n Motion, follow these steps:

1. Remove the device and other components from box.
2. Plug the *adapter cord* into the device socket and plug the *adapter* into an electrical outlet. All *energy level indicator lights* will turn on, and then go off.

9.1 Before your first treatment

Remove all jewellery before use.

You must decide which energy level is right for your skin, and the area you want to treat. The device has 5 levels of energy from 1 (the lowest) to 5 (the highest). For example, you might treat one area at energy level 1, and another at level 4. For information on energy levels, refer to section "8.2 Energy levels".

To choose the right energy level, test each treatment area separately (refer to the skin colour chart at the front of this manual).

1. Press the *control button* to turn on the device. The level 1 *energy level indicator light* will go on.
2. Place the *treatment surface* on your skin to deliver one pulse. If you do not experience abnormal discomfort press the *control button* again. This will increase the energy level. Deliver one pulse on a different spot at energy level 2.
3. Continue testing at increased energy levels. Use the energy levels indicated on the skin colour chart. When you find the highest energy level that is comfortable for that treatment area, you can start your treatment.

Use a lower energy level if you experience any discomfort.

Repeat this test for each part of the body you wish to treat.

You are now ready to use your Silk'n Motion

1. Before treating an area, the skin should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. **Do not wax, pluck or tweeze the hair out.**

Be sure to cover up moles, warts or other marks on the skin that should not be treated.

2. Press the *control button*. The device will turn on, the level 1 *energy indicator light* will be on, and a fan sound will start.
3. Approximately 1 second after pressing the *control button*, the *ready light* will turn on. The device is then ready to use at the lowest energy level.
4. **Set the energy level.** Press the *control button* repeatedly to increase the energy level until you have the level you want. The *energy level indicator lights* indicate the energy level setting (refer to section "8.2 Energy levels").
5. Decide if you wish to treat the area by **Gliding** or **Pulsing** (refer to section "8.3 Pulsing and gliding").
6. Apply the *treatment surface* to the skin. Make sure that the skin is spread evenly and smoothly. When the *treatment surface* is in full contact with the skin, the *ready light* will start blinking slowly.
7. The device will first determine the colour of your skin automatically and it will flash a pulse of light onto your skin. You will feel a mild sensation of warmth and tingling, and you may hear a popping sound.
8. If you are pulsing, remove the device from your skin. The device will immediately recharge for the next pulse. After a few seconds (depending on the energy level) the *ready light* will turn on again, and start blinking slowly when pressed against the skin. If you are pulsing, repeat the process. If you are gliding, simply continue.
9. Continue to pulse or glide until you have treated the area fully.

Stop use immediately, if your skin blisters or burns.

10. After treatment

When you have finished your treatment session:

1. Turn the device off by giving the *control button* an extended press. Remember the last energy level setting you used, because it will not be restored when you turn the device on again.
2. Remove the adapter from the electrical outlet.

11. Maintenance

It is important to clean your Silk'n Motion device after each treatment. To clean the device:

-  Unplug the device.
-  To clean the device including the *treatment surface*, wipe with a clean, dry cloth. Never use acetone or other abrasive fluids, because they will damage your device.
-  After cleaning, we recommend you store the device in its original box.
-  Never immerse your Silk'n Motion or any of its parts in water.

12. Troubleshooting

Do not attempt to open or repair your device . You may expose yourself to dangerous electrical components or pulsed light energy, and may cause serious harm to yourself. It will also void your warranty.

Problem	Checks
My device does not start.	Make sure the <i>adapter</i> is properly connected to the device and that it is plugged into an electrical outlet on the wall.
There is no light pulse when the device touches my skin.	Check that the device is connected and the <i>ready light</i> is on. Make sure that the <i>treatment surface</i> is evenly and firmly pressed to the skin and that you are moving the device across your skin.
One or more lights blink.	The device has malfunctioned.
The device has stopped working and is in standby mode.	Press the <i>control button</i> again to start the device.

For problems or queries, please contact Silk'n Customer Service.

13. Customer service

See our website for more information: www.silk'n.eu. This manual is also available as a PDF download from: www.silk'n.eu. Contact Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service number	Email
Belgium	0900-25006	serviceconsommateurbe@silk'n.eu
France	0891-655557	serviceconsommateurfr@silk'n.eu
Germany	089 51 23 44 23	kundenservice@inno-essentials.de
Italy		servizioclientiit@silk'n.eu
Netherlands	0900-2502217	klantenservicenl@silk'n.eu
Spain	807088200	servicioalcliente@silk'n.eu
UK	0906-2130009	customercareuk@silk'n.eu
Other countries	+31 (0)180-330 550	info@silk'n.com

14. Warranty

Refer to your Silk'n warranty card.

15. Disposal

Help protect the environment. Do not dispose with household waste and be sure to follow your country's laws for the disposal of electrical items.

16. Specifications

Model no.	H3220
Spot size	0.9cm x 3cm [2.7cm ²]
Technology	Home Pulsed Light
Maximum energy level	Max 4J/cm ²
Wavelength	475-1200nm
Operation & safety	Skin colour sensor ensures use on appropriate skin types.
Package size	220x220x135 mm
System weight	275 gm
Transport & storage between uses and storage condition	Temperature: -40°C to +70°C Relative humidity: 10% to 90% Atmospheric pressure: 500 to 1060 hPa
Operating conditions	Temperature: 10°C to 35°C Relative humidity: 30% to 75% Atmospheric pressure: 700 to 1060 hPa
Time from minimum storage temperature between uses to minimum operating temperature with ambient temperature of 20°C degrees	30 minutes
Time from maximum storage temperature between uses to maximum operating temperature with ambient temperature of 20°C degrees	30 minutes
Mode of operation	Continuous

Adapter models and rating	KSAPV0361200300D5 PMA36R-120
Input Output	100-240 Vac, 50-60 Hz, 0.9A 12 Vdc, 3.0A
Expected service life of the equipment	5 years

17. Signs and symbols

	Products marked with this symbol fulfil the requirements of the EU directives.
	Class II. The adapter is protected against electric shock by means of double insulation.
IP30	The device is protected against the ingress of objects larger than 2.5mm diameter, such as fingers, to protect against electric shock.
	Do not use in wet environments.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.

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